

Analysis Software* Quick Reference Guide

Run Details


Run Date: Tuesday, 1 January 2008

Run Number: Assay 1

Kit Batch Number: 053400401

Operator: ED A

Enter test information.

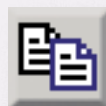

Paste Raw Data

Copy and paste data,
or enter values manually.

Raw Data

	1	2	3	4	5	6	7	8
A	1.970	0.103	3.012	1.126	1.199	0.071	0.075	0.166
B	0.050	0.623	1.416	0.355	0.364	0.102	0.400	0.253
C	0.053	0.199	2.460	0.112	0.131	0.660	0.795	1.432
D	0.052	0.193	1.379	0.068	0.051	0.062	0.062	1.989
E	0.053	0.079	0.906	0.064	0.073	0.052	0.195	3.001
F	0.065	0.167	1.108	0.056	0.432	0.063	0.063	0.248
G	0.074	0.097	0.265	0.050	0.047	0.075	0.079	0.262

Create and save format
(default or manual method).

Change Subject ID

Subject ID: A0090005

OK Cancel

ID Prefix: A009

OK Cancel

Enter subject names.
Calculate results.

Standards Results

Std	Conc	Mean	% CV	QC Result
S1	4.00	1.163	4.4	PASS
S2	1.00	0.360	1.8	PASS
S3	0.25	0.122	N/A	PASS
S4	0.00	0.060	N/A	PASS

Curve Details

Intercept: 0.9334 Slope: 0.8145

Correlation Coefficient: 1.00 PASS

ELISA Results
Valid ELISA test run.

View test Quality Control.


Subject Results

Subject ID	Nil	TB Ag	Mitogen	TB Ag- Nil	Mitogen- Nil	Result
A0090001	7.76	0.21	> 10	7.57	5.33	NEGATIVE
A0090002	0.09	1.89	5.19	1.80	5.10	POSITIVE
A0090003	0.09	0.47	> 10	0.38	> 10	POSITIVE
A0090004	0.09	0.45	5.02	0.36	4.93	POSITIVE
A0090005	0.09	0.15	3.00	0.06	2.91	NEGATIVE
A0090006	0.12	0.38	3.84	0.26	3.72	NEGATIVE
A0090007	0.14	0.19	0.66	0.05	0.52	NEGATIVE
A0090008	2.03	2.54	2.80	0.51	0.77	POSITIVE
A0090013	0.13	0.14	0.37	0.01	0.24	INDETERMINATE
A0090014	0.21	1.10	0.63	0.89	0.42	POSITIVE

View test results.



Print test results,
or save as PDF file.

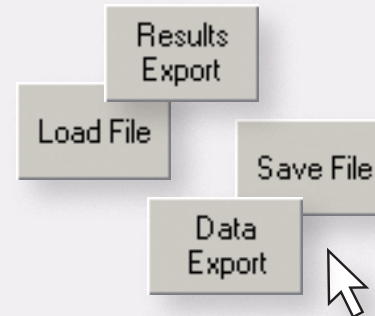

Further options.

For comprehensive instructions, please refer to the Analysis Software Instructional Guide

North America / South America

Email: customer.service@cellestis.com
Fax: +1 661 775 7479
Tel: +1 661 775 7480 (outside USA)
Toll free: 800 519 4627 (USA only)

Europe / Middle East / Africa

Email: europe@cellestis.com
Fax: +49 6151 428 59 110
Tel: +49 6151 428 59 0

Asia / Oceania

Email: quantiferon@cellestis.com
Fax: +61 3 9568 6623
Tel: +61 3 8527 3500